




vivo Surgical
Beyond Surgery

UNLOCK SUPERIOR SURGERIES

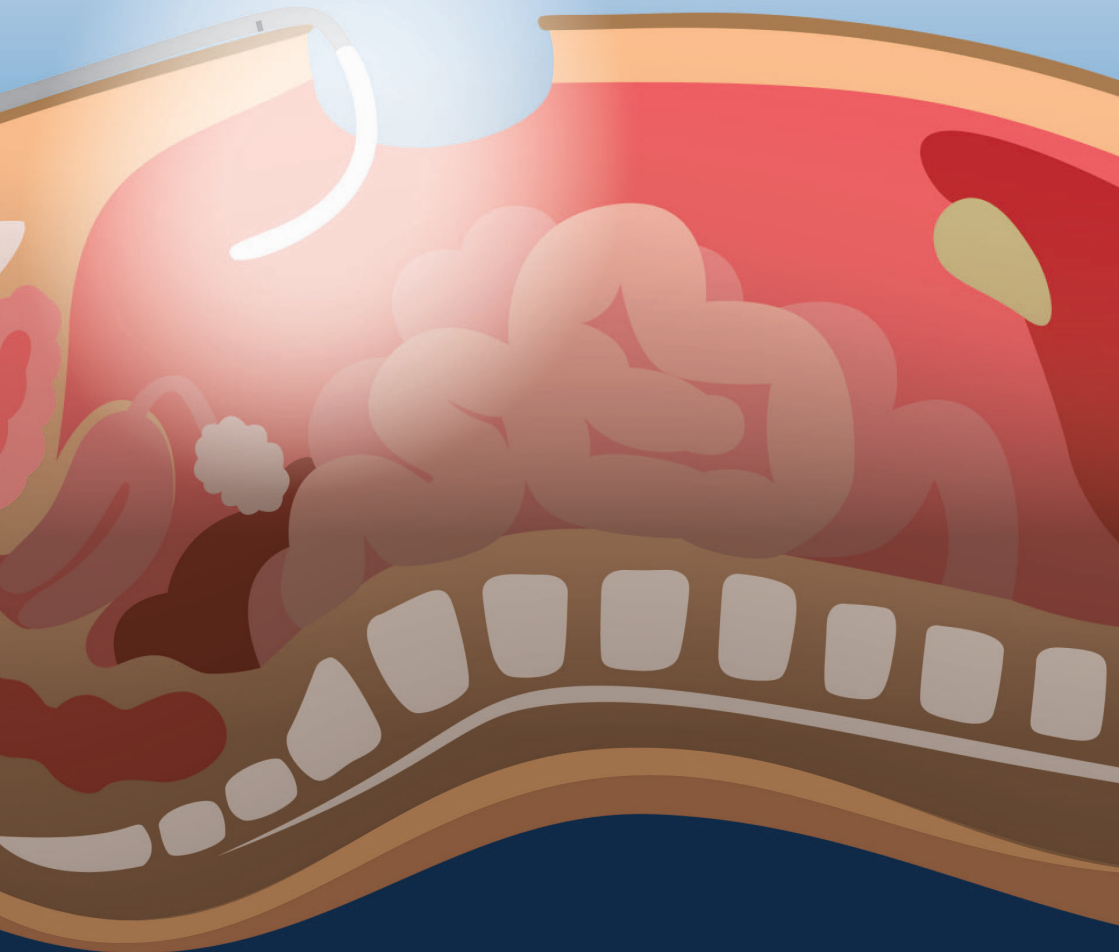

Deep Cavity In Vivo Surgical Lighting



KLARO™ is an innovative in vivo surgical lighting device that is far superior to other conventional lighting systems. The device optimally illuminates open surgical cavities, providing surgeons with a highly versatile and effective lighting solution especially suited for deep cavity work.



Adjustable Illumination

Choose from four lighting intensity levels to optimally illuminate deep surgical cavities



Cordless Convenience

Portable and lightweight with no wires to restrict movement



Long Battery Life

4-hour lifespan with integrated low-battery indicator (< 15 minutes left)



Safe Working Temperature

Maintains a safe working temperature under 38°C / 100.4°F within the surgical cavity at any lighting intensity



Sterile, Single Use

Supplied sterile hence eliminating risk of contamination



Wide-angle Lighting

When unbent: ~155° of illumination; when bent: up to ~340° of illumination



Ease of Use

Hassle-free set-up, operated with just a button push



Hands-free

The device clips onto surgical drapes or is placed on the surgical trolley during use



Small and Unobtrusive






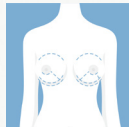


A light strip diameter of 4.6 mm ensures an easy fit in small, tight surgical cavities

Why Choose KLARO™

	KLARO™	Overhead Lighting	Headlamps	Lighted Retractors
1. Illuminates deep and narrow surgical cavities	V	-	-	-
2. Centres movement on the surgeon's immediate field	V	-	-	-
3. Reduces shadow effects within the surgical field	V	-	-	V
4. Reduces eye glare & fatigue for surgeons	V	-	-	-
5. No additional capital equipment necessary	V	-	V	-
6. Maintains safe working temperature	V	V	-	-

DID YOU KNOW?
 KLARO™ was invented by surgeons to solve key clinical needs in the surgical lighting space!

Applicable Surgeries

 Head and Neck	 General	 Orthopaedic	 Gynaecology
 Surgical Oncology	 Breast	 Cardiothoracic	 Oral & Maxillofacial

Transform **Any** Surgical Retractor into a **Lighted One**



KLARO™ can be easily secured onto any surgical retractor using our retractor loops.

Retractor loops are supplied with each KLARO™ device.

KLARO™ retractor loop instructional videos



For wide
handheld
retractors



For slim
handheld
retractors



For
self-retaining
retractors

Product Code	Description	Sales Unit
555-01-01	KLARO™ In Vivo Surgical Lighting Device (Sterile)	10 units per box

VIVOSurgical
Beyond Surgery

Address: 67 Ayer Rajah Crescent
#01-01/02
Singapore 139950

Enquiries: info@vivo-surgical.com
Sales: sales@vivo-surgical.com
Phone: +65 6677 0395

Follow us on: LinkedIn: VivoSurgical
Youtube: VivoSurgical
www.vivo-surgical.com